

RED LIGHT THERAPY BEDS

Red light therapy beds are full-body devices lined with thousands of medical-grade LEDs that emit specific red and near-infrared wavelengths (typically 630-940 nm). They deliver 360-degree coverage to stimulate cellular energy (ATP), reduce inflammation, relieve joint pain, and improve skin tone in 10-to-20-minute sessions.

**M4 | M4+**

COMPACT PLATFORM

**M5N | M5N+**

HIGH-OUTPUT POD

**M6N | M6N+**

PREMIUM FULL-BODY POD

**M7 | M7+**

STREAMLINED COMMERCIAL BED

WAVELENGTH CONTROL

Independent

PULSE SYSTEM

0-10,000 Hz

SESSION RANGE

10-20 min

COMMERCIAL COOLING

Integrated

Choose the Right Platform

Start with the room, electrical supply and expected appointment volume. Then choose the output level. Every family provides full-body red and near-infrared coverage, independent wavelength control, pulse modulation, Bluetooth music and commercial cooling.

FAMILY	BUYER POSITIONING	LEDs	IRRADIANCE	POWER	FOOTPRINT
M4 M4+	Compact and lighter	11,616 pcs 28,512 pcs	35-61 mW/cm ²	1,160 W 2,850 W	1935 x 855 x 990 mm
M5N M5N+	Wide, high-output pod	16,560 pcs 38,600 pcs	48-104 mW/cm ²	3,000 W 5,500 W	2275 x 1250 x 965 mm
M6N M6N+	Premium output and space	18,720 pcs 41,600 pcs	63-124 mW/cm ²	3,500 W 6,500 W	2275 x 1235 x 1125 mm
M7 M7+	Streamlined commercial bed	13,440 pcs 26,880 pcs	48-100 mW/cm ²	2,200 W 4,500 W	2145 x 1060 x 910 mm

SHARED COMMERCIAL FEATURES

CONTROL

Smart Controller 2.0 or optional Wireless Pad Controller 2.0.

LIGHT

Independent wavelength selection with continuous or pulsed operation.

COMFORT

Forced-air cooling, adjustable body fan and Bluetooth stereo music.

OPERATIONS

Preset protocols, manual mode, temperature monitoring and usage records.

SAFETY

Protective eyewear included; no ultraviolet light is used.

MOBILITY

Integrated wheels support installation, cleaning and service access.

M4 | M4+

A lighter full-body bed for space-conscious projects.

The M4 family combines a compact footprint, lower installed weight and a 110 V or 220 V option. It is a practical starting point for smaller wellness rooms, spas and private recovery spaces.



OPEN CONFIGURATION



CLOSED CONFIGURATION

WAVELENGTHS

633 / 660 / 850 / 940 nm

IRRADIANCE

35-61 mW/cm²

FOOTPRINT

1935 x 855 x 990 mm

WEIGHT

110 kg

Best fit: compact rooms, moderate appointment volume and projects that need a lower electrical load.

M4

LEDS

11,616 pcs

POWER

1,160 W

CURRENT

6 A at 220 V

M4+

LEDS

28,512 pcs

POWER


2,850 W

CURRENT


13 A at 220 V

M4 | M4+

PARAMETER	SPECIFICATION
Model	M4 M4+
Light source	Taiwan EPISTAR 0.2 W special LED
Wavelengths	633 / 660 / 850 / 940 nm
Irradiance	35-61 mW/cm ²
Wavelength control	Independent control
Pulse system	0-1,000 Hz
LED lifespan	Up to 100,000 hours
Control system	Smart Controller 2.0 / optional Wireless Pad Controller 2.0
Cooling	Self-circulating forced-air cooling
Bluetooth music	Yes
Input voltage	220 V or 110 V
Dimensions	1935 x 855 x 990 mm
Inner clearance	710 x 370 mm
Open height	1390 mm
Weight	110 kg
Body colors	White, Orange, Black
Eyewear / wheels	Included / Yes

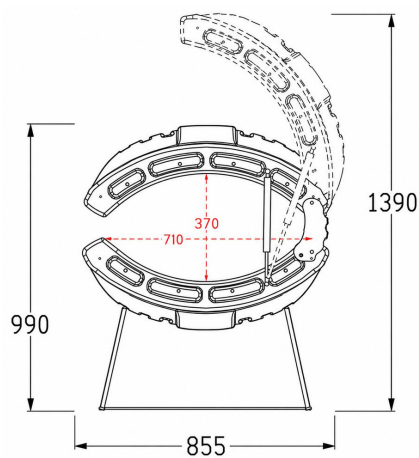

END VIEW - OPEN

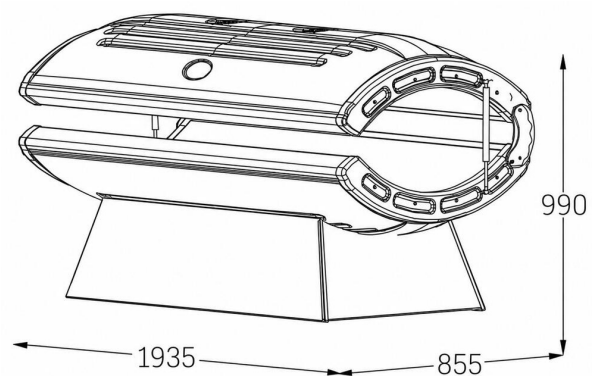

THREE-QUARTER VIEW - CLOSED

BUYER CHECK

Confirm delivery access, ceiling height, ventilation, voltage, circuit protection and the final plug or hard-wired connection before production.

M5N | M5N+

Five wavelengths with strong commercial output.

The M5N family is a wide, full-body pod for recovery rooms and wellness centers. It combines five independently controlled wavelengths, output from 63-124 mW/cm² and a commercial self-circulating cooling system.



OPEN CONFIGURATION



CLOSED CONFIGURATION

WAVELENGTHS

**633 / 660 / 810 / 850 /
940 nm**

IRRADIANCE

48-104 mW/cm²

FOOTPRINT

2275 x 1250 x 965 mm

WEIGHT

270 kg

Best fit: recovery rooms, spas and gyms planning regular full-body appointments.

M5N


LEDS	16,560 pcs
POWER	3,000 W
CURRENT	14 A at 220 V

M5N+


LEDS	38,600 pcs
POWER	5,500 W
CURRENT	25 A at 220 V

M5N | M5N+

PARAMETER	SPECIFICATION
Model	M5N M5N+
Light source	Taiwan EPISTAR 0.2 W special LED
Wavelengths	633 / 660 / 810 / 850 / 940 nm
Irradiance	48-104 mW/cm ²
Wavelength control	Independent control
Pulse system	0-1,000 Hz
LED lifespan	Up to 100,000 hours
Control system	Smart Controller 2.0 / optional Wireless Pad Controller 2.0
Cooling	Self-circulating forced-air cooling
Bluetooth music	Yes
Input voltage	220 V or 380 V
Dimensions	2275 x 1250 x 965 mm
Inner clearance	750 x 320 mm
Open height	1460 mm
Weight	270 kg
Body colors	Lunar White / Crystal Black
Eyewear / wheels	Included / Yes

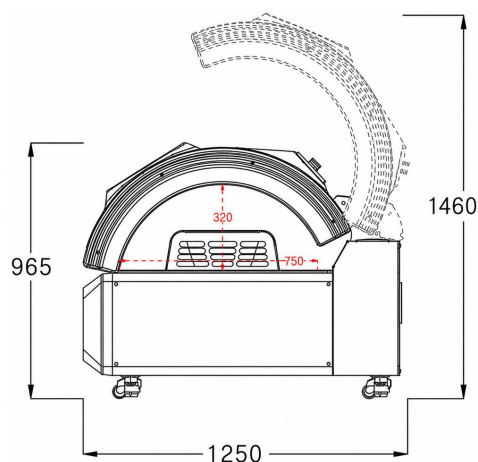

SIDE VIEW - OPEN

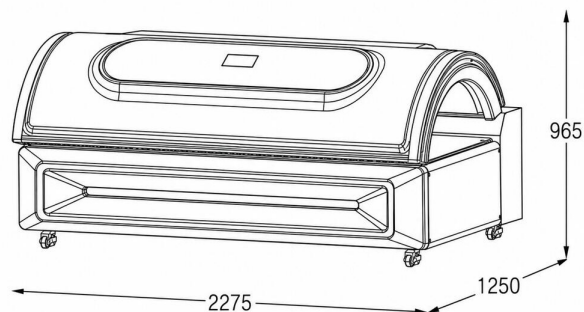

THREE-QUARTER VIEW - CLOSED

BUYER CHECK

Confirm delivery access, ceiling height, ventilation, voltage, circuit protection and the final plug or hard-wired connection before production.

M6N | M6N+

Maximum internal space and the highest rated output.

The M6N family adds side coverage and a larger inner clearance for premium red light therapy studios. It offers five wavelengths, irradiance from 48-104 mW/cm² and the highest LED count in the collection.



OPEN CONFIGURATION



CLOSED CONFIGURATION

WAVELENGTHS

**633 / 660 / 810 / 850 /
940 nm**

IRRADIANCE

63-124 mW/cm²

FOOTPRINT

2275 x 1235 x 1125 mm

WEIGHT

380 kg

Best fit: premium studios, high-specification projects and facilities expecting frequent use.

M6N


LEDS	18,720 pcs
POWER	3,500 W
CURRENT	16 A at 220 V

M6N+


LEDS	41,600 pcs
POWER	6,500 W
CURRENT	30 A at 220 V

M6N | M6N+

PARAMETER	SPECIFICATION
Model	M6N M6N+
Light source	Taiwan EPISTAR 0.2 W special LED
Wavelengths	633 / 660 / 810 / 850 / 940 nm
Irradiance	48-104 mW/cm ²
Wavelength control	Independent control
Pulse system	0-1,000 Hz
LED lifespan	Up to 100,000 hours
Control system	Smart Controller 2.0 / optional Wireless Pad Controller 2.0
Cooling	Self-circulating forced-air cooling
Bluetooth music	Yes
Input voltage	220 V or 380 V
Dimensions	2275 x 1235 x 1125 mm
Inner clearance	800 x 420 mm
Open height	1645 mm
Weight	380 kg
Body colors	Lunar White / Crystal Black
Eyewear / wheels	Included / Yes

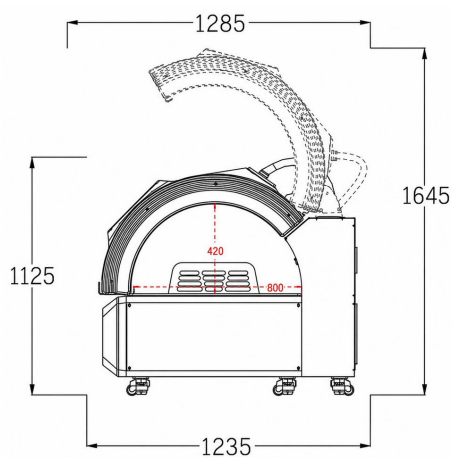

SIDE VIEW - OPEN

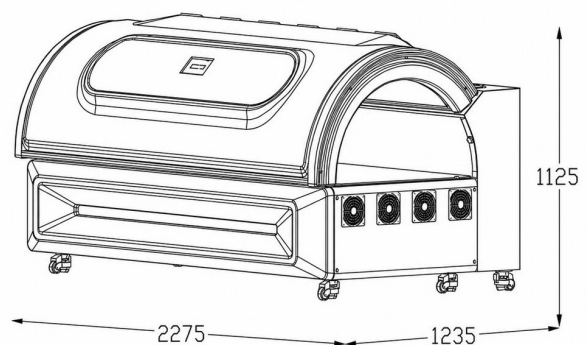

THREE-QUARTER VIEW - CLOSED

BUYER CHECK

Confirm delivery access, ceiling height, ventilation, voltage, circuit protection and the final plug or hard-wired connection before production.

M7 | M7+

Five-wavelength performance in a narrower profile.

The M7 family provides a streamlined commercial footprint with five independently controlled wavelengths. It is designed for spas, sports recovery rooms and wellness clinics that want a premium bed without the M6N footprint.



OPEN CONFIGURATION



CLOSED CONFIGURATION

WAVELENGTHS

**633 / 660 / 810 / 850 /
940 nm**

IRRADIANCE

48-100 mW/cm²

FOOTPRINT

2145 x 1060 x 910 mm

WEIGHT

270 kg

Best fit: spas, wellness clinics and sports recovery facilities balancing output and room size.

M7


LEDS	13,440 pcs
POWER	2,200 W
CURRENT	12 A at 220 V

M7+


LEDS	26,880 pcs
POWER	4,500 W
CURRENT	21 A at 220 V

M7 | M7+ / Planning

PARAMETER	SPECIFICATION
Model	M7 M7+
Light source	Taiwan EPISTAR 0.2 W special LED
Wavelengths	633 / 660 / 810 / 850 / 940 nm
Irradiance	48-100 mW/cm ²
Wavelength control	Independent control
Pulse system	0-1,000 Hz
LED lifespan	Up to 100,000 hours
Control system	Smart Controller 2.0 / optional Wireless Pad Controller 2.0
Cooling	Self-circulating forced-air cooling
Bluetooth music	Yes
Input voltage	220 V or 380 V
Dimensions	2145 x 1060 x 910 mm
Inner clearance	810 x 370 mm
Open height	1435 mm
Weight	270 kg
Body colors	Lunar White / Crystal Black
Eyewear / wheels	Included / Yes

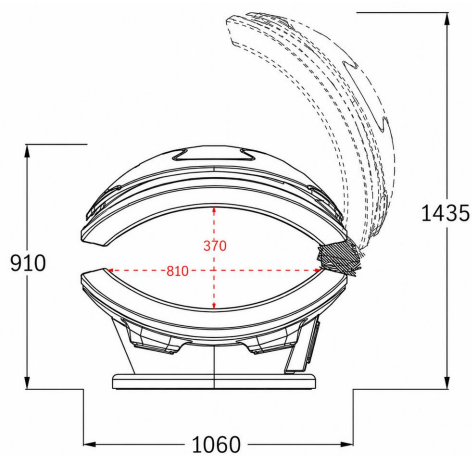

END VIEW - OPEN

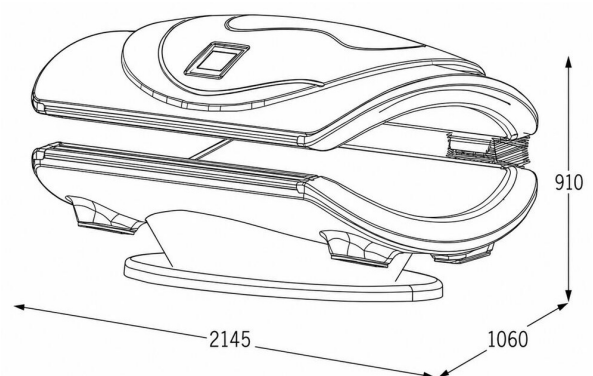

THREE-QUARTER VIEW - CLOSED

BUYER CHECK

Confirm delivery access, ceiling height, ventilation, voltage, circuit protection and the final plug or hard-wired connection before production.

Specify the Project Before the Model

A successful installation starts with practical project information. The following five checks prevent most specification changes after production has started.

01

ROOM ENVELOPE

Measure the treatment room, doorway, corridor, lift and canopy opening height. Leave working clearance for cleaning and service.

02

ELECTRICAL SERVICE

Confirm voltage, phase, dedicated circuit capacity, breaker, cable size and grounding with a qualified local electrician.

03

APPOINTMENT VOLUME

Estimate realistic daily sessions. Higher output may suit frequent use, but the room, cooling and electrical load must support it.

04

OPERATING WORKFLOW

Decide whether staff will use preset protocols, manual settings, external tablet control or centralized operating rules.

05

MARKET CONFIGURATION

Confirm body color, controller, logo, voltage, electrical connection, packaging and certification requirements in the PI.

TYPICAL SESSION WORKFLOW**PREPARE**

Expose the intended area and wear eyewear.

SELECT

Choose the approved preset or manual protocol.

RUN

Typical sessions are 10-20 minutes.

CLEAN

Clean the acrylic after every client.

COOL

Allow the fan-delay cycle to finish.

BUILD THE ROOM AROUND THE BUYER.

Send the project details first. We will help match the platform, electrical configuration and control option to the intended facility.

INCLUDE THESE DETAILS WITH YOUR RFQ

- Preferred model or family
- Destination country and voltage
- Available room dimensions
- Expected daily appointments
- Controller and OEM requirements
- Delivery access and installation plan

Shenzhen Medlight Technology Co., Ltd

Room 306, Bld D, Dingfeng Scientific Park, Dalang Town, Longhua District, Shenzhen, China 518109

Web www.medledlight.com

Cell +86-13554796007

Email sales@medledlight.com